


## **Press Release: Service Support Ends for discontinued Rancimat**

Date: 8th<sup>th</sup> July 2010

Service support for the 679 Rancimat ended at the end of 2009 due to the model being discontinued.

In addition, the following accessories have been discontinued

- Glass jointed reaction vessels
- Glass jointed measuring vessels
- Glass conductivity cells

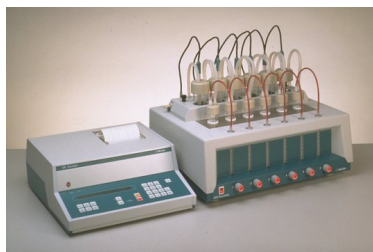

*The discontinued 679 Rancimat*

### **EXCLUSIVE OFFER!**

Because we can no longer support the 679, until the end of 2010 Metrohm will offer all existing 679 users a discount of 25% of the list price of the 743 Rancimat.

### **743 RANCIMAT**

The 743 Rancimat provides a type of measuring arrangement with disposable reaction vessels and evaluation software with comprehensive database functions which allows the handling of large amounts of samples and data. The 743 Rancimat has two separately temperature controlled heating blocks, each with up to four reaction vessels. Up to four instruments can be run by a single PC which means it is possible to run 32 samples at 8 different temperatures in parallel. A GLP set allows qualification of the measuring system. As all measuring vessels, reaction vessels and gas tubing are now consumables; this eliminates contamination and reduces cost when replacing glassware.



2.743.0014 743 Rancimat for oils and fats including software and accessories

For a quotation and further information please contact us on 01928 579 600 or contact us on [enquiry@metrohm.co.uk](mailto:enquiry@metrohm.co.uk)